

## DNA Marker V

### PRODUCT FORMATION

Components	Specification	Cat. No.
DNA Marker V	250 $\mu$ l	MD1012
DNA Marker V	250 $\mu$ l $\times$ 5	MD1112

### STORAGE

The product can be stored at normal temperature (0-30 $^{\circ}$ C) for more than three years. If the product is not used for a long period of time, please store at -20 $^{\circ}$ C to prevent the evaporation of water.

### TECHNICAL SUPPORT


TEL: 400-0099-857

E-MAIL: [technical@simgen.cn](mailto:technical@simgen.cn)

### INTRODUCTION

DNA Marker V is composed of 7 individual DNA fragments, presenting 5.5k, 3.5k, 2.5k, 1.5k, 1k, 500, 250 bp sharp bands respectively. DNA Marker V contains 1 $\times$ Loading Buffer, users can apply 5 - 10  $\mu$ l in agarose gel electrophoresis directly.

The red and yellow tracking dye in DNA Marker V will not weakened the DNA bands under UV light, better than bromophenol blue and xylene cyanol FF.

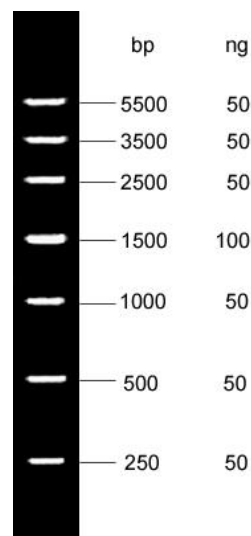

### PRECAUTION

1. Clear bands can be obtained by applying 5  $\mu$ l DNA Marker when the lane width is less than 5 mm. If the lane is wider, loading volume of DNA Marker should be increased appropriately.
2. For DNA electrophoresis, agarose purity is of great significance to DNA band definition. Therefore, agarose with good quality should be used and gel concentration of 1% ~ 2% is recommended.
3. During agarose electrophoresis, the concentration of agarose is closely associated with the separation of DNA fragments. High agarose concentration is ideal for the separation of the short DNA fragments. While low agarose concentration is ideal to separate the long DNA fragments.